

REPORT DOCUMENTATION PAGE			Form Approved OMB No. 0704-0188	
Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503.				
1. AGENCY USE ONLY (Leave blank)		2. REPORT DATE 22.Apr.02	3. REPORT TYPE AND DATES COVERED MAJOR REPORT	
4. TITLE AND SUBTITLE NANOMETER SCALE ANTIBODY PATTERNING FOR DIRECTED CELL IMMOBILIZATION AND STIMULATION			5. FUNDING NUMBERS	
6. AUTHOR(S) CAPT ORTH REID N				
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) CORNELL UNIVERSITY			8. PERFORMING ORGANIZATION REPORT NUMBER CI02-58	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES) THE DEPARTMENT OF THE AIR FORCE AFIT/CIA, BLDG 125 2950 P STREET WPAFB OH 45433			10. SPONSORING/MONITORING AGENCY REPORT NUMBER	
11. SUPPLEMENTARY NOTES				
12a. DISTRIBUTION AVAILABILITY STATEMENT Unlimited distribution In Accordance With AFI 35-205/AFIT Sup 1			12b. DISTRIBUTION CODE	
13. ABSTRACT (Maximum 200 words)				
20020523 157				
14. SUBJECT TERMS			15. NUMBER OF PAGES 5	
			16. PRICE CODE	
17. SECURITY CLASSIFICATION OF REPORT	18. SECURITY CLASSIFICATION OF THIS PAGE	19. SECURITY CLASSIFICATION OF ABSTRACT	20. LIMITATION OF ABSTRACT	

#026328

SECURITY AND POLICY REVIEW WORKSHEET- Request for Public Release Clearance (In Accordance With AFI 35-101, Chapter 15) (Do not use this worksheet to request clearance of software or web pages.)			
SUBMITTING ORGANIZATION (Office Symbol)		AFIT/CIGG	
A. DOCUMENT TYPE (Documents must be complete including all figures, charts, photographs and text.)			
<input type="checkbox"/> ABSTRACT	<input type="checkbox"/> JOURNAL ARTICLE	<input type="checkbox"/> THESIS	<input type="checkbox"/> DISSERTATION
<input type="checkbox"/> BROCHURES	<input type="checkbox"/> NEWS RELEASE	<input type="checkbox"/> SPEECH	<input checked="" type="checkbox"/> OTHER
<input type="checkbox"/> CD-ROM	<input type="checkbox"/> PHOTO WITH CAPTIONS	<input type="checkbox"/> SUCCESS STORY	Conference Paper
<input type="checkbox"/> DISPLAY/EXHIBIT	<input type="checkbox"/> POSTER SESSION	<input type="checkbox"/> TECHNICAL PAPER	
<input type="checkbox"/> FACT SHEET	<input type="checkbox"/> PRESENTATION WITH TEXT	<input type="checkbox"/> TECHNICAL REPORT	
B. TITLE OF DOCUMENT Nanometer-Scale Antibody Patterning for Directed Esophageal Cell Immobilization and Stimulation			NO. PAGES 5
C. AUTHOR(S) NAME AND DUTY TITLE Capt Reid Orth, AFIT Ph.D. Candidate			OFFICE SYMBOL AFIT/CIGG
D. FORUM (Public release clearance is not required for material presented in a closed meeting and which will not be made available to the general public on the Internet or in any print or electronic media)			
DOCUMENT WILL BE <input type="checkbox"/> PRESENTED ORALLY	NAME OF CONFERENCE	LOCATION	DATE
DOCUMENT WILL BE <input checked="" type="checkbox"/> DISTRIBUTED IN PRINT	NAME OF PUBLICATION IEEE-EMAS Special Topics Conference Proceedings		SUBMITTAL DEADLINE ASAP
HAVE RELATED DOCUMENTS BEEN PREVIOUSLY CLEARED FOR PUBLIC RELEASE <input type="checkbox"/> NO <input checked="" type="checkbox"/> YES AFIT CASE NUMBER			
NOTE: If document will be released to a medium outside the Department of Defense, the following disclaimer must be added: The views expressed in this article are those of the author and do not reflect the official policy or position of the United States Air Force, Department of Defense, or the U.S. Government.			
FOR AFIT/PA USE ONLY		SUBMITTED: Air Force Institute of Technology Publications AFIT/PA Building 642 Room 2010 2950 P Street Wright-Patterson AFB OH 45433-7765	
		FOR MORE INFORMATION Phone: (937) 255-9354 DSN: 785-9354 FAX: (937) 255-2135 DSN: 785-2135 E-MAIL: afit.pa@afit.edu WEB: http://www.afit.mil/newsrm01.htm	

E. NATIONAL SECURITY AND TECHNOLOGY ISSUES

CHECK YES OR NO (Provide a brief explanation for all positive responses in the SYNOPSIS SECTION below.)

- ☐ YES ☒ NO DOES THE MATERIAL HAVE THE POTENTIAL TO BECOME AN ITEM OF NATIONAL OR INTERNATIONAL INTEREST.
- ☐ YES ☒ NO DOES THE MATERIAL AFFECT NATIONAL SECURITY POLICY OR FOREIGN RELATIONS.
- ☐ YES ☒ NO DOES THE MATERIAL CONCERN SUBJECTS OF POTENTIAL CONTROVERSY AMONG DOD COMPONENTS OR WITH OTHER FEDERAL AGENCIES.
- ☐ YES ☒ NO DOES THE MATERIAL CONTAIN TECHNICAL DATA DEVELOPED UNDER CONTRACT OR INDEPENDENTLY DEVELOPED AND CONTROLLED BY THE INTERNATIONAL TRAFFIC IN ARMS REGULATIONS (ITAR) THAT MAY BE MILITARILY CRITICAL AND SUBJECT TO LIMITED DISTRIBUTION, BUT ON WHICH A DISTRIBUTION DETERMINATION HAS NOT BEEN MADE.
- ☐ YES ☒ NO DOES THE MATERIAL CONTAIN INFORMATION ON NEW WEAPONS OR WEAPON SYSTEMS, SIGNIFICANT MODIFICATIONS OR IMPROVEMENTS TO EXISTING WEAPONS OR WEAPON SYSTEMS, EQUIPMENT OR TECHNIQUES.
- ☐ YES ☒ NO DOES THE MATERIAL CONTAIN INFORMATION ON NATIONAL COMMAND AUTHORITIES; COMMAND, CONTROL, COMMUNICATIONS, COMPUTERS AND INTELLIGENCE; INFORMATION WARFARE; OR COMPUTER SECURITY.
- ☐ YES ☒ NO DOES THE MATERIAL CONTAIN INFORMATION ON MILITARY ACTIVITIES OR APPLICATIONS IN SPACE, NUCLEAR WEAPONS, INCLUDING WEAPON-EFFECTS RESEARCH; CHEMICAL AND BIOLOGICAL WARFARE ISSUES; BIOLOGICAL AND TOXIN RESEARCH, HIGH-ENERGY LASERS AND PARTICLE BEAM TECHNOLOGY; ARMS CONTROL TREATY

F. SYNOPSIS (Provide a brief description of the system, process, or technology, including its state of development and whether the application is military, commercial, or dual-use). EXAMPLE: This is new concept of applying current development for high-powered lasers, specifically for military applications.

Biological Materials (protein, eosinophil cells, antibodies) have been patterned through a stencil. Immune system cells are reacted to this surface. This is a primary research tool, for commercial applications.

G. REFERENCES (Check all that apply)

- ☒ ALL REFERENCES ARE UNCLASSIFIED, UNLIMITED, AND ARE AVAILABLE TO THE PUBLIC
- ☐ REFERENCE # ARE SUBJECT TO DISTRIBUTION LIMITATION. NO LIMITED INFORMATION FROM THESE REFERENCES IS INCLUDED IN THE DOCUMENT
- ☐ NO REFERENCES ARE INCLUDED IN THIS DOCUMENT

RELEASE AUTHORITY AND MANAGEMENT RESPONSIBILITY FOR JOINT EFFORTS

IS THIS INFORMATION CO-AUTHORED WITH SOMEONE FROM AN ORGANIZATION OUTSIDE OF AFIT.

☒ YES ☐ NO

DOES YOUR ORGANIZATION HAVE THE AUTHORITY TO RELEASE THE INFORMATION

☒ YES ☐ NO

IF WORK IS SPONSORED BY ANOTHER ORGANIZATION, DO YOU HAVE THE AUTHORITY TO RELEASE ALL THE INFORMATION

☒ YES ☐ NO

IF NOT, WHAT ORGANIZATION SHOULD COORDINATE ON THE REQUEST

ORGANIZATION

PHONE

NOTE: If the organization is a DOD organization, AFIT/PA will coordinate release approval with that organization's PA office

1. AUTHOR/TECHNICAL REVIEW CERTIFICATION	
AUTHOR	
I certify the information contained in the attached document is technically accurate and does not disclose classified, sensitive, or militarily critical technology and does not violate proprietary rights or copyright restrictions. All security and technology issues listed above have been considered and any applicable security classification guides have been reviewed.	
SIGNATURE <i>Reid N. Orth</i>	DATE 22 Mar 02
AUTHOR(S) NAME (Print) REID. N. ORTH	TELEPHONE 607-255-6286
OFFICE SYMBOL/ORGANIZATION AFIT/CIGG	E-MAIL ADDRESS rno2@cornell.edu
TECHNICAL REVIEWER CERTIFICATION (Faculty Advisor)	
I certify the information contained in the attached document is technically accurate and does not disclose classified, sensitive, or militarily critical technology and does not violate proprietary rights or copyright restrictions. All security and technology issues listed above have been considered and any applicable security classification guides have been reviewed.	
SIGNATURE	DATE
TECHNICAL REVIEWER'S NAME (Print)	TELEPHONE
OFFICE SYMBOL/ORGANIZATION	E-MAIL ADDRESS
DEPARTMENT HEAD (or equivalent)/PROGRAM MANAGER CERTIFICATION	
The attached material submitted for public release has been properly staffed and reviewed by this division/organization. I support the recommendation of the Author and Technical Reviewer for public release.	
SIGNATURE <i>Melissa L. Flattery</i>	DATE 26 MAR 02
NAME (Print) MELISSA L. FLATTERY	TELEPHONE 5-3152x3026
OFFICE SYMBOL/ORGANIZATION AFIT/CIGG	
AFIT PUBLIC AFFAIRS	
THE ATTACHED MATERIAL IS NOT CLEARED FOR PUBLIC RELEASE	
SIGNATURE OF SECURITY AND POLICY REVIEW OFFICER <i>Pat Wagner</i>	DATE 28 Mar 02
J. REMARKS (please reference the specific section)	

Nanometer-Scale Antibody Patterning for Directed Cell Immobilization and Stimulation

Reid N. Orth, Lori Kwan, Ismail Hafez, Jun Kameoka, Manfred Lindau, and Harold G. Craighead

School of Applied and Engineering Physics, Cornell University, Ithaca, New York 14853

Antibodies have been patterned at nanoscale resolution for the precise immobilization and stimulation of immunological cells. We demonstrate that an antigen, bovine serum albumin (BSA), can be patterned on silicon using a photolithographically patterned polymer lift-off technique. The nanoscale pattern is realized as the polymer is mechanically peeled away in one contiguous piece in solution. Anti-BSA antibodies bound specifically to the BSA create a pattern of antibody F_c segments that provide the stimulus for eosinophils immobilization and degranulation. These patterns, ranging from sub 600 nm to 67 μm , provide a spectrum of stimuli for the 10-14 μm eosinophil cells. This method provides a new technique for capturing cells from solution, analyzing cellular biochemical cascades events such as degranulation, and studying cellular morphological changes in response to a finite, nanoscale antigenic stimulus.

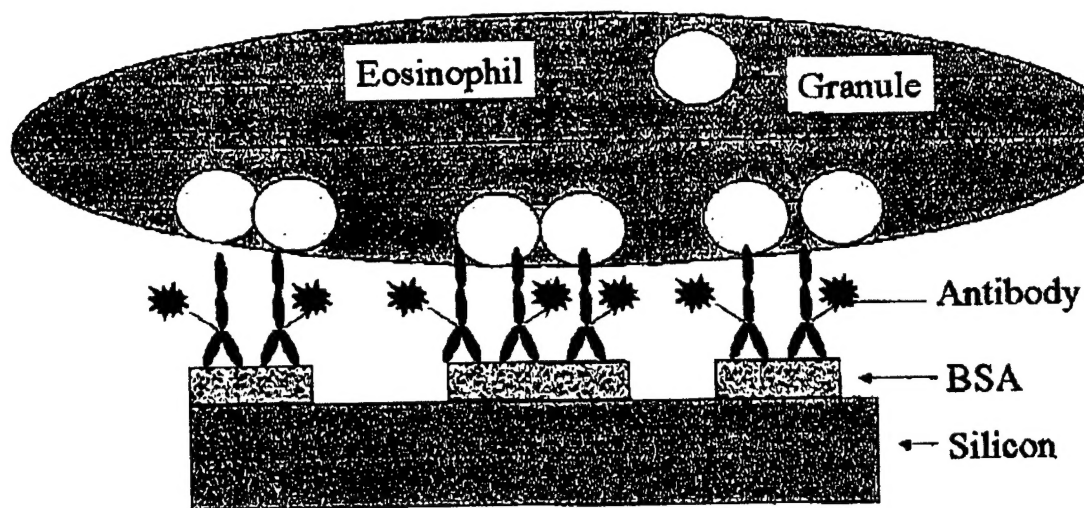


Figure 1. Schematic of eosinophil cell binding to Fc tail of Alexa 594 conjugated anti-BSA antibody on silicon substrate.

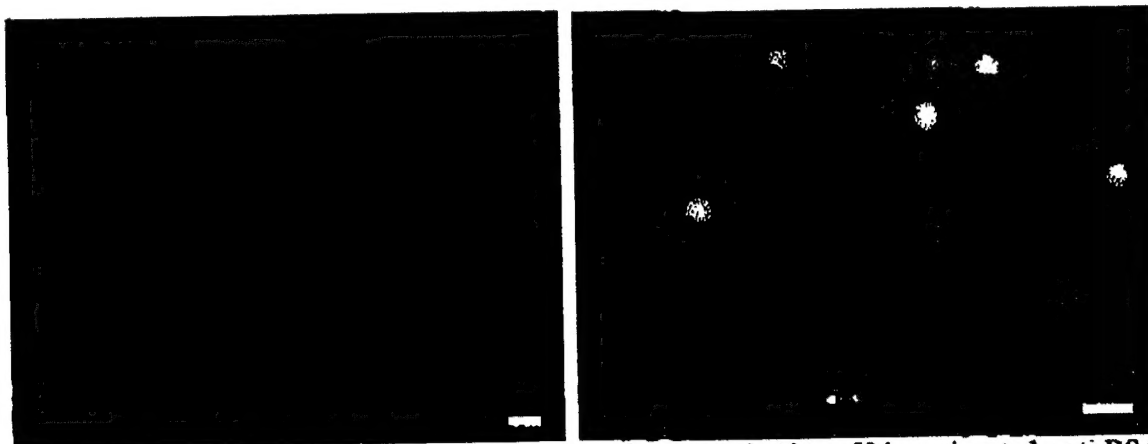


Figure 2A and 2B. (A) Epifluorescent image of patterned Alexa-594 conjugated anti-BSA antibodies specifically bound to patterned BSA on silicon substrate. Scale bar is 70 μm and the largest square is 67 μm . (B) Eosinophil cells (green) immobilized to the patterned antibody matrix. Scale bar is 19 μm and the period between dots is 9.6 μm .

NANOMETER-SCALE ANTIBODY PATTERNING FOR DIRECTED CELL IMMOBILIZATION AND STIMULATION

Reid N. Orth, Lori Kwan, Jun Kameoka, Manfred Lindau, and Harold G. Craighead

School of Applied and Engineering Physics, Cornell University, Ithaca, New York 14853

Abstract - Antibodies (Ab) are patterned at nanoscale precision for the precise immobilization and stimulation of immune cells. We demonstrate that the antigen bovine serum albumin (BSA) can be patterned on silicon using a photolithographically patterned polymer lift-off technique. The nanoscale pattern is realized as the polymer is mechanically peeled away in one contiguous piece in aqueous solution. Anti-BSA Ab are bound specifically to BSA to create a pattern of oriented Ab that provides a surface for eosinophil immobilization and degranulation. The patterns ranged from $0.36 \mu\text{m}^2$ to $4,489 \mu\text{m}^2$, appropriate dimensions for the 10-14 μm diameter eosinophil cells. This method provides a new technique for immobilizing cells onto nano- and micrometer scale patterns for analyzing cellular biochemical cascade events such as degranulation and studying cellular morphological changes in response to defined nanoscale antigenic stimulus.

Keywords - Immunology, nanofabrication, biomaterials, surface modification.

1. INTRODUCTION

Eosinophilic granulocytes (eosinophils) play an important role in the host immune defense against parasitic invasion [1]. These exocytotic cells synthesize cytotoxic proteins and store these proteins in secretory granules, where they are ready to be released onto the surface of a parasite [2]. Eosinophils contribute to the pathophysiology of bronchial asthma [3] with activation of allergen-specific IgG1 and IgG3 induced through the Fc portion of the immunoglobulin molecules (Fc γ R1) [4]. Eosinophil cell Fc and complement receptor numbers vary dramatically between healthy patients and patients with eosinophilia [5]. Eosinophils are stimulated and degranulate in response to zymosan (yeast cell walls) [6], IgG-coated Sephadex beads [7], IgA coated beads [8], and lipid mediators leukotriene B₄ (LTB₄) and platelet activating factor (PAF) [9]. Whole cell patch-clamp technique has been used to study the discrete electrical impedance changes during compound exocytosis and cumulative fusion [10]. Cellular adhesion through the B₂ integrin is an important step in eosinophil activation and accumulation as demonstrated by eosinophil binding to IgG coated sepharose beads and by binding inhibition with membrane antibodies against CD18 and CD11b--the B₂ integrin ligands [11].

Micron-scale patterning of biomolecules is an active area of research as the microelectronics technology merges with biology [12]. Patterning biomaterials on the micro- and nanometer scales allow a more focused method for cell

stimulation than bath application of the biomaterial. Several different methods have been used to micropatterned biomaterials and chemicals on solid substrates for cell interrogation. Photopatterning has been used to spatially distribute biomolecules, such as enzymes, antibodies, and nucleic acids, for the development of biochips on silicon, glass, and plastic substrates [12]. A printing technique called microcontact printing (μ CP) process uses a poly(dimethyl siloxane) (PDMS) elastomeric stamp to pattern a wide array of and biomaterials [14],[15],[16]. Microfabrication techniques were used to micropattern bovine serum albumin and horseradish peroxidase [17]. μ CP has been used to pattern cell adhesion proteins to immobilize and direct neuronal and astroglial cell growth on glass substrates [18],[19]. Deep plasma etching and photoplastics were used to create PDMS stamps to pattern different cell suspensions to specific locations of a tissue culture substrate [19].

Many different cell types have been exposed to patterned biomaterials on solid substrates. *E.coli* O157:H7 cells were captured from a solution using PDMS stamped anti-*E.coli* O157:H7 for a diffraction grating biosensor [20]. Neuronal and glia cells were also patterned using light-assisted functionalized photoresists [22]. A silicon micromachined flow-through chamber has been designed for the 'entrapment' of chick embryo spinal cord neurons as a model system for biological neural networks [23]. Differentiated B104 neuroblastoma cells were micropatterned on four substrates to determine the preferred support substrate [24].

This paper presents a method for patterning a surface with antibodies that serve as a stimulus for eosinophil immobilization, activation and degranulation. A patterned BSA:anti-BSA IgG complex serves as a model stimulus for eosinophils. Vapor deposited Parylene, di-para-xylylene, is conformally deposited over the silicon substrate and used as a pinhole-free barrier between the solution and the substrate. Conventional photolithography and reactive ion etching (RIE) are used to pattern the polymer [25]. After the BSA incubation, the samples are incubated with Alexa-594-conjugated anti-BSA IgG solution. The immunospecificity of the antigen allows Alexa-594-conjugated anti-BSA IgG molecules to bind specifically to the patterned substrate. The samples are subsequently immersed in buffer solution and the Parylene layer is mechanically removed. Micro- and nanometer square patterns are formed using a polymer lift-off method [26] (Fig. 1). Eosinophils are pipetted onto

the patterned surface (Fig. 2). Time series epifluorescence microscopy occurs immediately after cell addition to detect real time cell immobilization, morphological changes, and degranulation.

II. METHODOLOGY

1) *Silicon Wafer Preparation and Parylene Deposition:* 3-inch <1-0-0> N/phos type wafers (Silicon Qwest, Int'l, Santa Clara, CA) are cleaned in base and acid baths to remove surface contaminants. The wafers are baked at 1100°C for 50 minutes and annealed for 10 minutes in a wet oxide process to grow a 500-nm thermal oxide layer in the silicon substrate. Fig. 1 details the fabrication steps used for the parylene lift-off technique. A pinhole-free conformal layer of Parylene is deposited onto 3-inch silicon wafers using the PDS-2010 Labcoater 2 Parylene deposition system (Specialty Coating Systems, Indianapolis, IN). The polymer thickness is dependent upon the amount of evaporated polymer. 1.5 g of Parylene C dimer is used to deposit a 1- μ m thick Parylene film on five 3-inch silicon wafers.

2) *Photolithography:* 1.5 μ m of OCG_OiR 897-12i photoresist (Shipley, Marlboro, MA) is applied to the Parylene-coated silicon wafers. The samples are pre-baked for 1 minute at 90°C and exposed using standard photolithographic techniques with a 10X stepper (Fig. 1a). After development, the exposed portions of the Parylene film are postbaked for 90 seconds at 115°C and are subjected to an oxygen-based RIE step using the PI asma Therm 72 with an RF power density at 0.255 W/cm² (Fig. 1b). After etching, the samples are dipped into a beaker of acetone to remove residual photoresist, rinsed with isopropyl alcohol, and washed in deionized water. The samples are then dried with a nitrogen gas stream.

3) *Silanization of Silicon:* 3-aminopropyltriethoxysilane (3-APTS, Sigma-Aldrich, Milwaukee, WI) solution is prepared in a 50-mL amber bottle using 0.5-mL of 3-APTS and 24.0 mL of acetone in a nitrogen purged glovebox to create a 2% silane solution. The silanization step began by cleaning 1 cm² silicon chips in a Harrick Plasma Cleaner/Sterilizer PDC 3-G for 1 minute. The chips are removed and placed in 100°C Milli-Q filtered water for 30 minutes. The silicon chips are nitrogen dried then quickly inserted into the bottled silane solution and incubated in a closed container for 30 minutes. The chips are removed, washed in acetone for 5 minutes, immersed in isopropyl alcohol and deionized water, and baked on a hotplate at 70°C.

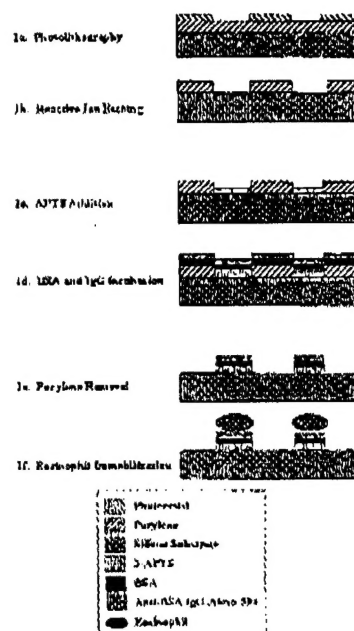


Fig. 1. Process flow schematic of the fabrication steps. (a) Patterning of 1.5 μ m of OCG_OiR 897-12i photoresist using optical lithography. (b) Reactive ion etching of 1.0 μ m layer of Parylene C dimer and subsequent removal of the top photoresist layer. (c) Application of 3-aminopropyltriethoxysilane layer to a plasma cleaned silicon substrate. (d) Application of 100 μ g/ml HSA and 50 μ g/ml anti-BSA IgG. (e) Mechanical parylene removal with tweezers, resulting in a surface with patterned Alexa 594-BSA-anti-BSA complex. (f) Immobilization of eosinophils to the patterns IgG surface.

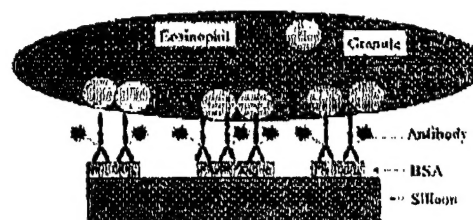


Fig. 2. Schematic of eosinophil cell binding to Fc tail of Alexa 594 conjugated anti-BSA antibody on silicon substrate.

4) *BSA and Anti-BSA IgG Preparation:* Bovine Serum Albumin (BSA) (Sigma-Aldrich, Milwaukee, WI) is reconstituted to 10 mg/ml and used to create 100, 10, and 1 μ g/ml dilutions. A 30 μ L drop of 2 mM BSA solution is placed on the Parylene-patterned substrate for 60 minutes, as illustrated in Fig. 1c. Polyclonal, mouse anti-BSA IgG molecules (Sigma-Aldrich, Milwaukee, WI) and stained with NHS-Alexa 594 dye (Molecular Probes, Eugene, OR). The IgG stock solution is diluted to 100 μ g/ml, 7.4 pH in phosphate buffered saline (PBS). Samples are incubated in

35-mm plastic Petri dishes (Fisher Chemicals, Pittsburgh, PA). Alexa 594-conjugated anti-BSA IgG is applied onto the pattern and incubated for 30 minutes. After the incubation, the sample remains immersed in aqueous solution while being transferred to a second Milli-Q water beaker in a 35-mm Petri dish. The Parylene is removed mechanically by peeling it off the substrate with tweezers in solution. The polymer film is shed easily in one contiguous piece from the substrate. The resulting sample contains patterned antigen-Ab complexes as illustrated in Fig. 1d.

5) Eosinophil Cell Preparation and Immobilization: Fresh blood is drawn from the jugular veins of horses, eosinophils are isolated, and are purified over discontinuous Percoll gradients as previously described [27]. Purified eosinophils are suspended in Medium 199 (Sigma-Aldrich, Milwaukee, WI) containing 4 mM glutamine, 4.2 mM NaHCO₃ and penicillin/streptomycin (Sigma-Aldrich, Milwaukee, WI) (pH 7.2-7.3) stored at room temperature and used within 24 hours. The cell suspension is incubated for 45 minutes in LysoTracker Green (Molecular Probes, Eugene, OR). The samples are placed in the bottom indentation of 35mm petri dish (Mat Tek, Ashland, MA) containing 2 ml of PBS buffer. 50 μ l of 1.9×10^6 cells/ml solution is pipetted directly onto the patterned silicon chips. Time series epifluorescence imaging begins immediately after cell application using an Olympus AX 70 upright microscope with a 60X water immersion objective and Omega Optical filter sets (Brattleboro, VT). Alexa 594 dye is observed with a 510-590-nm excitation/590-nm emission filter set and LysoTracker Green is observed with a 450-490-nm excitation/520-nm emission filter set. Images are captured using a Spot CCD camera (Diagnostic Instruments, Inc., Sterling Heights, MI). Cells are fixed with 3.7% paraformaldehyde (Sigma, Milwaukee, WI) for long-term storage.

III. RESULTS

In this report, we describe the micro- and nanometer scale arrays of patterned Ab for cell stimulation, test the patterned surface antigenicity with functional antibodies, develop a method for immobilizing eosinophil cells to the pattern, analyze the cells' preferential binding to the patterned surface, and detect eosinophil degranulation events after cell incubation on the pattern.

Fig. 3a illustrates Alexa 594-conjugated anti-BSA IgG squares with sub 600-nm widths on the lower right to 67- μ m squares on the upper left) on a silicon substrate using the Parylene lift-off technique. This matrix of patterns demonstrates the wide range of sizes attainable with this patterning technique. 500 nm resolution is approximately the lowest photopatterning threshold attainable using the 10X stepper. A statistical analysis showed that the average relative aerial density of fluorescently labeled antigen to the background is at least 150 times greater in the exposed regions than in the lift-off regions [28].

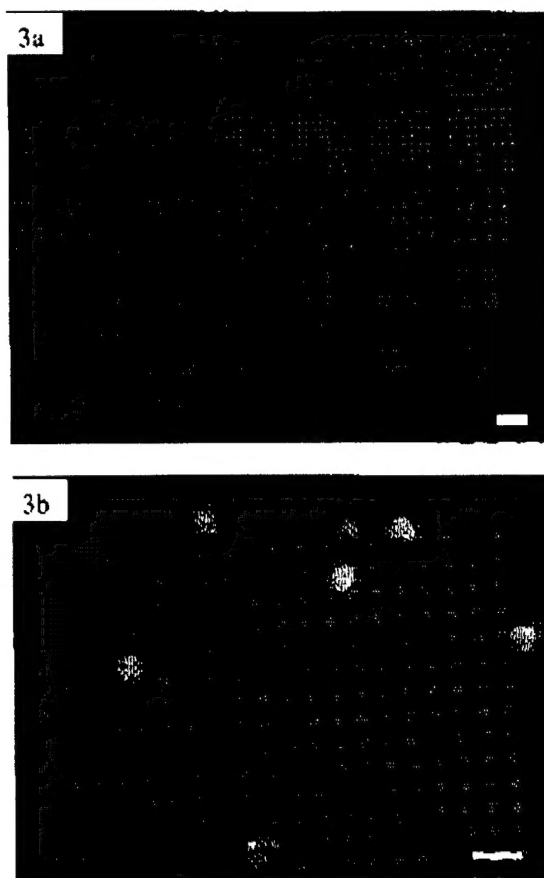


Fig. 3a and 3b. (a) Epifluorescent image of patterned Alexa-594 conjugated anti-BSA antibodies specifically bound to patterned BSA on silicon substrate. Scale bar is 70 μ m and the largest square is 67 μ m. (b) Eosinophil cells (green) immobilized to the patterned antibody matrix. Scale bar is 19 μ m and the period between dots is 9.6 μ m.

Fig. 3b illustrates the immobilization of horse eosinophils to the BSA-anti-BSA IgG complexes. Horse eosinophils were used since they have large granules that are easy to visualize during microscopy. This binding is similar to the experiments using immunoglobulin coated Sepharose beads detailed in the introduction. This method provides information that Sepharose beads do not offer. First, this method can be used to determine the spatial distribution and size of patterned Fc stimuli required for eosinophil activation and degranulation. Second, this technique allows eosinophil migration to be analyzed over the planar patterned surface.

Fig. 4a shows the initial state and Fig. 4b shows the same region after 5 minutes. These images demonstrate clear morphological changes are occurring on the substrate. Consequently, the circled regions in this image show possible fusion and/or degranulation events may be

occurring. Further experimentation with patch clamp amperometry will be required to confirm that these changes are taking place.

IV. DISCUSSION

This technique may be useful to determine if fusion events and degranulation events may be occurring on the patterned substrate using fluorescence imaging. The signal observed in unpatterned regions is comparable to the signal observed on a blank, new silicon wafer surface. Thus, the unpatterned regions show negligible binding. Uniformity of the final patterned biomaterial relies upon an optimized photoresist thickness of 1.5 μm and Parylene thickness 1.0 μm , photolithography precision, optimized RIE duration to prevent under- and overetching, and sufficient biomaterial incubation time. Application of 1 μm of Parylene and 1.5 μm of photoresist provided optimal conditions for polymer removal and resolution.

Anti-BSA Ab are used to confirm the specificity of the binding onto the patterned BSA and provide a surface rich in Fc fragments onto which eosinophils can bind and be stimulated. The patterned antigen did not spread from the confined pattern regions and were still isolated after a month of storage in aqueous solution.

The Parylene lift-off technique offers a rapid and precise way to create supported micro- and nanometer-scale patterns. These patterns can be used to capture other biomaterials from solution and integrated into biosensors, bioMEMS, and biological assay systems.

There are several advantages of this technique. The Parylene film removal can be performed at any step during the processing. Therefore, multiple reagents can be added to the initial patterned surface prior to Parylene removal. This would allow subsequent reagents to be added at high concentrations and maximal binding without concern for nonspecific binding on unpatterned surface areas. Parylene is biologically compatible polymer, provides a conformal coating with low permeability, and can be removed with a one-step mechanical lift-off. This technique permits the sample to remain submerged in solution so that the functional molecules are not denatured when dried. Parylene does not have a permanent bond with the substrate, thus allowing easy removal in one piece. The conformal film of Parylene is pinhole-free, so no unwanted patterning occurs in unexposed regions.


V. CONCLUSION

We have demonstrated a method for precise patterning antigen at the micro- and nanometer scale for cell stimulation. The antigenicity is confirmed by binding fluorescent antibodies onto the patterned surface. Eosinophil cells are effectively immobilized on the patterned surface. The eosinophils demonstrated preferential binding to the patterned surface. Time course

imaging provides a means to study spatial-temporal changes associated with stimulation and degranulation.

This technique offers a new versatile tool to pattern antigen and other biomaterial onto solid substrates with feature sizes below 600 nm.

4a

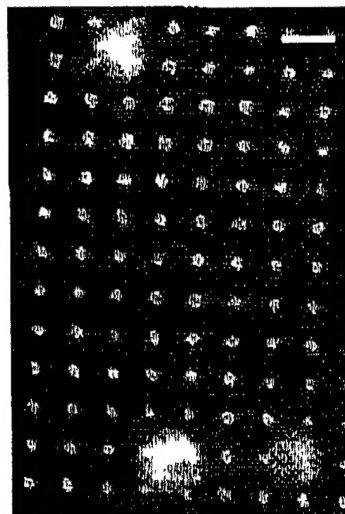

4b

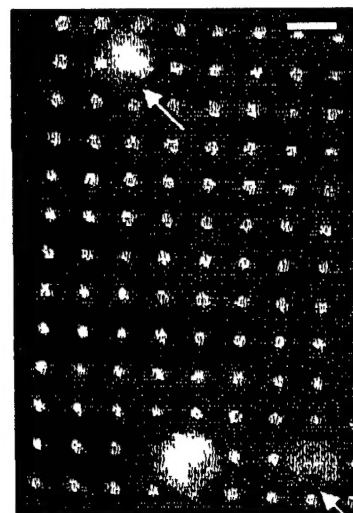


Fig. 4a and 4b. Epifluorescence images of eosinophils interacting with the patterned substrate. (a) Eosinophil interaction with the substrate 1 minute after application. (b) Eosinophil interaction with the substrate 6 minutes after application. Morphological changes can be observed in the cells with white circles. The cell with the top arrow appears to lose a green granule at its base. The cell with the lower arrow appears it may have degranulated from its center onto the pattern. Patch clamp amperometry will be required to confirm that these changes are taking place. Scale bar is 15 μm .

ACKNOWLEDGEMENT

We acknowledge the support of DARPA, NSF support through the Nanobiotechnology Center, and the resources of the Cornell Nanofabrication Facility. We would like to thank Dr. Ismail Hafez for his assistance preparing lipids and guidance on configuring a lipid extrusion system. We would also like to thank Mathieu Foquet for his designing the microscope set-up for this experiment and guidance with imaging. The views expressed in this article are those of the authors and do not reflect the official policy or position of the United States Air Force, Department of Defense, or the U.S. Government.

REFERENCES

- [1] J. Klein and V. Horejff. *Immunology*. Oxford: Blackwell Science Ltd., 1998, p. 52.
- [2] S. Sceppek, R. Mougel, and M. Lindau. "Compound exocytosis and cumulative degranulation by eosinophils and their role in parasite killing." *Parasitology Today*, vol. 10, p. 276, 1994.
- [3] J.H. Butterfield, K.L. Leiterman, and G.J. Gleich. "Eosinophil-associated diseases." in *M. M. Samter's Immunological Diseases*, vol. 1, M.M. Frank, K.F. Austen, H.N. Claman, and L.R. Umhoe, Eds. Boston: Little, Brown and Company, 1994, pp. 501-527.
- [4] M. Kaneko, M.C. Swanson, G.J. Gleich, and H. Kita. "Allergen-specific IgG1 and IgG3 through Fc gamma RII induce eosinophil degranulation." *J Clin Invest*, vol. 95, pp. 2813-21, June 1995.
- [5] I. Winqvist, T. Olofsson, I. Olsson, A.M. Persson, and T. Hallberg. "Altered density, metabolism and surface receptors of eosinophils in eosinophilia." *Immunology*, vol. 47, pp. 531-539, November 1982.
- [6] G. Archer and J. Hirsch. "Motion picture studies on degranulation of horse eosinophils during phagocytosis." *J. Experimental Medicine*, vol. 118, pp. 287-313, 1963.
- [7] I. Winqvist, T. Olofsson, and I. Olsson. "Mechanisms for adherence of eosinophils to an antibody-coated surface." *Scand J Haematol*, vol. 35, pp. 88-95, July 1985.
- [8] R.I. Abu-Ghazaleh, T. Fujisawa, J. Mestecky, R.A. Kyle, and G.J. Gleich. "IgA-induced eosinophil degranulation." *J Immunol*, vol. 142, pp. 2393-2400, April 1989.
- [9] C. Kroegel, T. Yukawa, G. Dent, P. Chanez, K.F. Chung, and P.J. Barnes. "Platelet-activating factor induces eosinophil peroxidase release from purified human eosinophils." *Immunology*, vol. 64, pp. 569-561, July 1988.
- [10] M. Lindau, J. Hartmann, and S. Sceppek. "Three distinct fusion processes during eosinophil degranulation." *Annals of the New York Academy of Sciences*, vol. 710, pp. 232, March 1994.
- [11] S. Kaneko, S. Horie, M. Kato, G.J. Gleich, and H. Kita. "A crucial role for b2 Integrin in the activation of eosinophils by IgG." *J Immunology*, vol. 155, pp. 2631-2541, 1995.
- [12] D.W. Branch, R.C. Wheeler, G.J. Brewer, and D.E. Leckband. "Long-term maintenance of patterns of hippocampal pyramidal cells on substrates of polyethylene glycol and microstamped polylysine." *IEEE Transactions on Biomedical Engineering*, vol. 47, pp. 290-300, 2000.
- [13] N. Dontha, W.B. Nowall, and W.G. Kuhr. "Generation of biotin/avidin/enzyme nanostructures with maskless photolithography." *Anal Chem*, vol. 69, pp. 2619-2625, July 1997.
- [14] H.A. Biebuyck, N.B. Larsen, E. Delamar, and B. Michel. "Lithography beyond light: Microcontact printing with monolayer resists." *IBM J. Res. Develop.*, vol. 41, pp. 159-170, 1997.
- [15] Y. Xia, and G.M. Whitesides. "Soft lithography." *Angewandte Chemie*, vol. 37, pp. 550-575, 1998.
- [16] P. Yang, T. Deng, D. Zhao, P. Feng, D. Pinc, D.F. Chmelka, G.M. Whitesides, and G.D. Stucky. "Hierarchically ordered oxides." *Science*, vol. 282, pp. 2244-2246, 1998.
- [17] S. Britland, E. Perez-Arnaud, P. Clark, B. McGinn, P. Connolly, and G. Moures. "Micropatterning proteins and synthetic peptides on solid supports: A novel application for microelectronics fabrication technology." *Biotechnol Prog*, vol. 8, pp. 155-160, March 1992.
- [18] C.D. Jones, R.C. Davis, L. Kani, H.G. Craighead, M. Isaacson, J.N. Turner, and W. Shain. "Patterned protein layers on solid substrates by thin stamp microneedle." *Langmuir*, vol. 14, pp. 741-744, 1998.
- [19] L. Kam, W. Shuin, J.N. Turner, and R. Bizios. "Correlation of astroglial cell function on micro-patterned surfaces with specific geometric parameters - hepatocytes and 3T3 fibroblast." *Biomaterials*, vol. 20, pp.2343-2350, 1999.
- [20] A. Folch, A. Ayon, O. Hurtado, M.A. Schmidt, and M. Tancr. "Molding of deep polydimethylsiloxane microstructures for microfluidics and biological." *J Biomech Eng*, vol. 121, pp. 28-34, 1999.
- [21] P.M. St. John, R. Divis, N. Cady, J. Czajka, C.A. Bart, and H.G. Craighead. "Diffraction-based cell detection using a microcontact printed antibody grating." *Anal Chem*, vol. 70, pp. 1108-1111, 1998.
- [22] D.V. Nicolau, T. Tsuchi, H. Taniguchi, H. Tanigawa, and S. Yoshikawa. "Patterning neuronal and glia cells on light-assisted functionalized photoreists." *Biosens Bioelectron*, vol. 14, p. 317, 1999.
- [23] S. Martinovic, M. Gove, M. Tedesco, D. Murgess, and M. Granatino. "A simple microfluidic system for patterning populations of neurons on silicon micromachined substrates." *J Neurosci Methods*, vol. 87, pp. 35-44, 1999.
- [24] J.M. Corey, A.L. Brunette, M.S. Chen, J.A. Weyhenmeyer, G.J. Brewer, and B.C. Wheeler. "Differentiated B104 neuroblastoma cells are a high-resolution assay for micropatterned substrates." *J Neurosci Methods*, vol. 75, pp. 91-98, 1997.
- [25] B. Ilic and H.G. Craighead. "Topographical patterning of chemically sensitive biological materials using a polymer-based dry lift off." *Biomedical Microdevices*, vol. 2, pp. 317-322, 2000.
- [26] R. Orth, T.G. Clark, and H.G. Craighead. "Avidin-biotin micropatterning methods for biosensor applications." In review.
- [27] S. Sceppek and M. Lindau. "Focal exocytosis by eosinophils - compound exocytosis and cumulative fusion." *EMBO J*, vol. 12, pp. 1811-1817, May 1993.
- [28] R. Orth, B. Ilic, T.G. Clark, and H.G. Craighead. "Patterning micro- and nanometer scale lipid bilayers using a polymer-based lift-off." Unpublished.